

FONA – Every dental solution you need

Through decades of experience and deep understanding of the dental profession, we deliver complete, reliable and accessible solutions. Regardless of country or specialisation, every dentist can find a fitting solution in FONA today. As our knowledge

and expertise grow, so does our range. Our international team is constantly on the move, meeting partners and dentists, so we can understand their needs and develop products that support and enhance their daily work.

Advanced dental technology

DIGITAL DENTISTRY



CHAIRSIDE CAD/CAM



EXTRAORAL IMAGING



INTRAORAL IMAGING

TRADITIONAL DENTISTRY



TREATMENT CENTERS



INSTRUMENTS AND LASER



HYGIENE

Headquarters

FONA Dental s.r.o.
Ševčenkova 34
85101 Bratislava
Slovak Republic

info@fonadental.com
www.fonadental.com

DISTRIBUTOR

Campeon

TAILOR-MADE TO FIT YOUR NEEDS





Campeon

TAILOR-MADE TO FIT YOUR NEEDS

Campeon is an elegant and highly advanced treatment center, created to meet the expectations of even the most demanding dental professionals. With emphasis on high-end finishing and latest technology, it is truly a champion in its field.

01

PATIENT COMFORT


Smooth chair movement and design guarantee comfort of the patient during treatment. Body contoured backrest and headrest ensure proper positioning with minimal stress to the spine.



Swingable armrests offer grip support when needed and allow easy reach of the treatment area.



Available with a 3D headrest for optimal patient comfort and convenient access during demanding treatments.

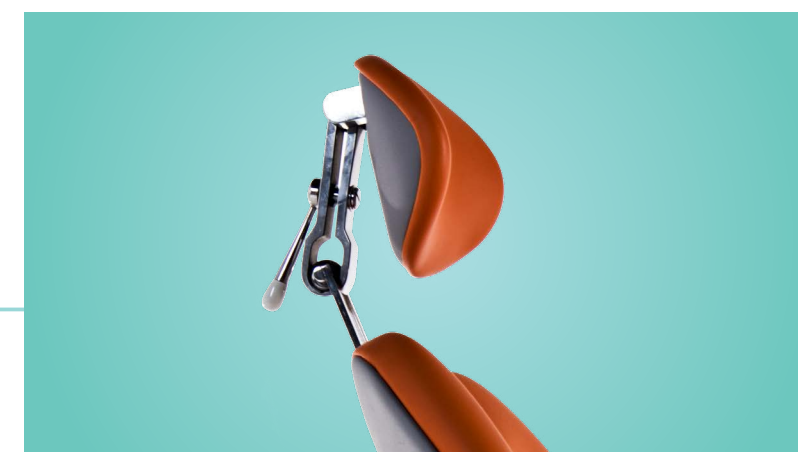

02

SUPERIOR ERGONOMY




Designed with dentist comfort in mind, Campeon is equipped with pantographic arms, swingable water unit and upholstered backrest among other work-enhancing features. It is also available in a version for left-handed users.

The pantographic assistant arm can be equipped with 3 to 5 instruments, advanced control panel and can be easily adjusted to desired treatment position.



Two foot control models are available with different functionalities to reflect your individual needs and preferred work position.

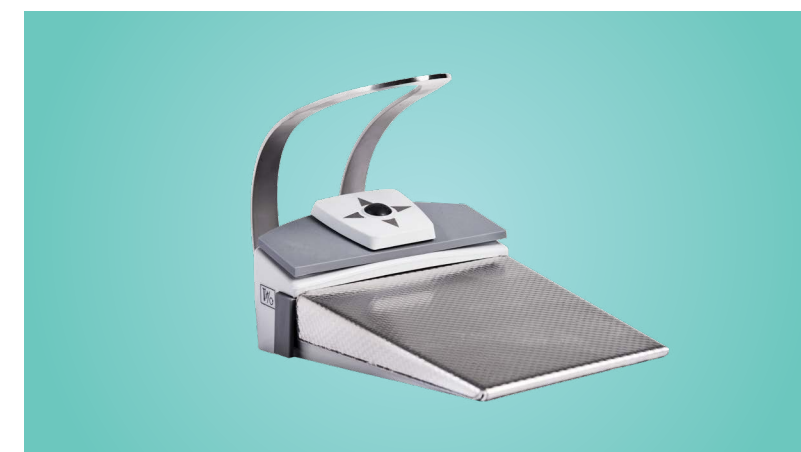

03

FEATURE VARIETY



Campeon is available in whip arm and hanging hose dentist element versions. The smart control panel allows to store up to 8 settings for each rotary instrument or scaler among a wide range of functions.

The dentist element can be equipped with up to 5 instruments, including endo motor, curing light, scaler and high speed handpieces with light.



Safety instruments interlock function, regular turbine maintenance reminder and programmable individual instrument cooling saves time and effort.




04

EXCEPTIONAL DESIGN

Campeon is made from quality materials, using latest production methods, resulting in a beautiful design with high-end finish and emphasis on every detail.

Choose from 11 colors and seamless or soft upholstery to perfectly fit the look of your dental office.

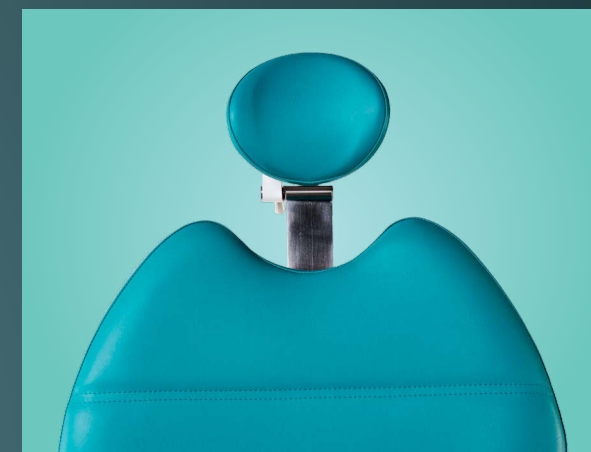

Variety of comfortable dentist and assistant stools is available to suit your preferred work style.



05

SMART HYGIENE AND SAFETY




Campeon offers a choice of disinfection systems and a rounded design to avoid cross-contamination. Multiple safety switches ensure protection for you and your patient during treatment.

There are three different disinfection systems to choose from. Semi-automatic or automatic for water lines and automatic for suction lines.



To ensure your safety there are anti-collision switches located on the backrest, assistant element and base of the chair.

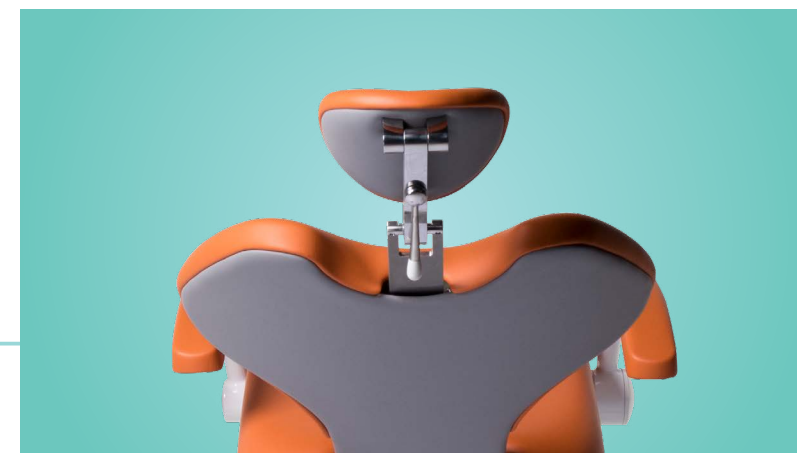

FEATURES AND OPTIONS



LIGHT

Sensor-operated LED Light with intensity up to 30,000 Lux, 3 axes of movement and removable silicone handles.



ASSISTANT ELEMENTS

Advance or Pro assistant element can be equipped with 3 to 5 instruments.



MONITOR

Optional 17" or 22" Neovo dental monitor for a diagnosis easily shareable with the patient.



TRAY

Tray on top or below the dentist element depending on version selected.



BOWL ROTATION

Ceramic cuspidor bowl swivelling towards the patient for comfortable access.



SMART DISPLAY

Embedded display allows full control over the setting of the instrument used.

TECHNICAL SPECIFICATION

MAX. POWER INPUT AT 230V/50HZ	400 VA
INPUT PRESSURE OF AIR	from 0,45 to 0,8 MPa
INPUT PRESSURE OF WATER	from 0,3 to 0,6 MPa
DENTAL UNIT WEIGHT (WITHOUT THE CHAIR) ADVANCE H	100 kg + max. 20 kg acc. to version
DENTAL UNIT WEIGHT (WITHOUT THE CHAIR) ADVANCE W	105 kg + max. 20 kg acc. to version
MAX. LOADING CAPACITY OF THE TRAY TABLE	1,5 kg

AVAILABLE COLORS



GRN4	TYR1	LBU1	BLU6	DBU1	GRY1
PLU2	VI02	VI03	ORA3	YEL3	